


TO DETERMINE **MELTING POINTS**



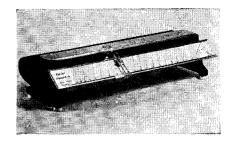

KOFLER HOT BENCH

Developed primarily to determine quickly the melting points and eutectic temperatures of organic substances, the Kofler Hot Bench also provides a rapid method to observe the presence and behavior of water crystal, the process of decomposition and to judge volatility.

The metal body is $1\frac{1}{2}$ " wide x $14\frac{1}{2}$ " long and is electrically heated at one end to provide an almost linear temperature drop. The temperature ranges from 50° to 260°C. The scale has temperature lines marked at every 2° and is read by a runner with pointer and tab. Test substances with known melting points are provided to gauge or set the instrument.

To determine a melting point, the substance is laid directly on the chromeplated surface. With pure substances, there is a sharp dividing line between the solid and fluid stage. The runner is moved until the pointer rests on the dividing line and the temperature is read from the scale. The margin of error is usually not more than \pm 1°C.

Available for prompt delivery FOR DETAILED INFORMATION WRITE FOR CATALOG AB

National Representatives: WM. J. HACKER & CO., Inc. 82 Beaver St., New York 5, N. Y.

CONTENTS, continued

Precision of Simple Flame Photometer
Interlaboratory Study on Determination of Acetyl in Cellulose Acetate 400
MICROCHEMISTRY
Determination of Microgram Quantities of Sulfide Sulfur
Spectrographic Determination of Beryllium in Urine and Air
Photometric Determination of Microquantities of Arsenic
NOTES ON ANALYTICAL PROCEDURES
Determination of Free Carbon in Atmospheric Dust
Separation of 2,4-Dinitrophenylhydrazones of Some Substituted Benzaldehydes by Chromatographic Adsorption
Mass Spectrometric Determination of Oxygen in Water Samples 414 Israel Dostrovsky and F. S. Klein
Rapid Determination of Organic Acids in Cured Tobacco
Technique for Noting Combustion in Oxygen Bombs
Rapid Volume and Concentration Calculations in Industrial Air Analysis 417 P. W. Mullen
Ascorbic Acid in Analytical Chemistry
Chicago Chemical Conference
Symposium on Modern Methods of Analytical Chemistry 425
Corrections
CRYSTALLOGRAPHIC DATA Ethylenedinitramine
SCIENTIFIC COMMUNICATION
BOOK REVIEW
THE ANALYST'S CALENDAR
AIDS FOR THE ANALYST Wave Spreader for Precision Determination of Half-Wave Potential with Sargent Model XII Polarograph
Thermostated H-Cell for Polarography
Adapters for Introducing Hot Coils into Distillation Flasks for Constancy of Ebullition
THE ANALYST'S COLUMN
INSTRUMENTATION
NEW PRODUCTS
MANUFACTURERS' LITERATURE

The American Chemical Society assumes no responsibility for the statements and opinions advanced by contributors to its publications. Views expressed in the editorials are those of the editors and do not necessarily represent the official position of the American Chemical Society.

Published by the American Chemical Society, from 20th and Northampton Sts., Easton, Pa. Executive Offices, Editorial Headquarters, and Circulation Department, 1155 Sixteenth St., N.W., Washington 6, D. C. Advertising Office: 332 West 42nd St., New York 18, N. Y. Branch Editorial Offices: Chicago 4, III., Room 819, 25 East Jackson, Blvd., Houston 2, Tex., 623 West Bldg., London, Bush House, Aldwych, New York 17, N. Y., 60 East 42nd St., San Francisco 2, Calif., Room 454, 750 Market St. Entered as second-class matter at the Post Office at Easton, Pa., U.S.A., February 18, 1948, under the Act of March 3, 1879, as 12 times a year monthly on the 15th. Acceptance for mailing at special rate of postage provided for in Section 1103, Act of Oct. 3, 1917, authorized July 13, 1918. Annual subscriptions: members \$3.50, nonmembers \$4.00. Postage to countries not in the Pan-American Union \$1.20, Canadian postage \$0.35. Combination rate for ANALYTICAL CHEMISTRY and Industrial and Engineering Chemistry: members \$7.00, nonmembers \$8.00. Postage for the combination: countries not in the Pan-American Union \$3.90, Canadian postage \$1.10. Single copies: current issues, members \$0.40, nonmembers \$0.50, back numbers, prices on request, special rates to members. Claims for missing numbers will not be allowed if received more than 60 days from date of issue. No claims allowed from subscribers in central Europe, Asia, or the Pacific islands other than Hawaii, or because of failure to notifythe Circulation Department of a change of address, or because copy is "missing from files."

CHANGE OF ADDRESS: Notify the American Chemical Society at 1155 Sixteenth St., N. W., Washington 6, D. C., of any change of address. Such notification should include both old and ne